

Experience advanced

Negative Pressure Wound Therapy

RELIABLE | EFFECTIVE | SEAMLESS | ECONOMICAL



Switch to
Smart NPWT Ecosystem


NoWoundTM
Smart Healing

by MEDVITAL VENTURES

NoWound™ Advanced NPWT

Accelerates wound healing through the following mechanisms:

- A Wound contraction**
NoWound aNPWT creates a vacuum effect that gently pulls the wound edges together, reducing the wound size and speeding up the healing process.
- B Microdeformation**
NoWound aNPWT strains the healing tissues, promoting cell division, growth factor production, and angiogenesis.
- C Fluid removal**
NoWound aNPWT removes excess fluid from the wound, reducing edema and improving blood flow through vasomotor mediators.
- D Reduced bacterial colonization**
NoWound aNPWT may provide a barrier that shields the wound from environmental contaminants.

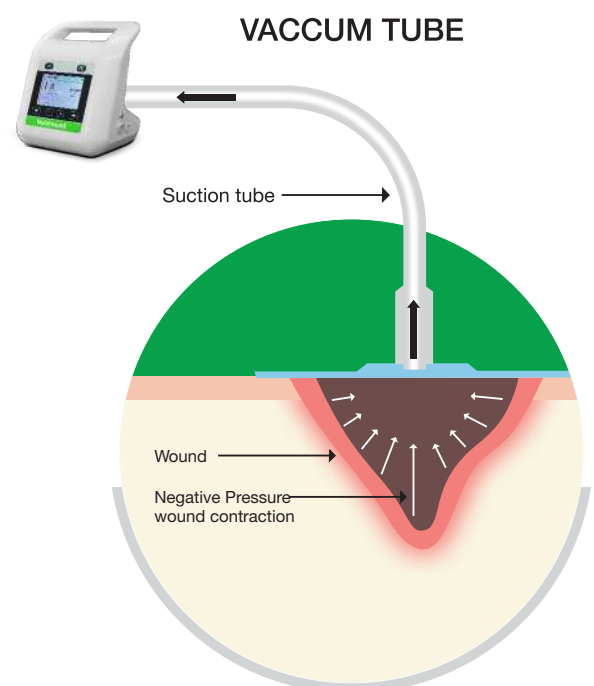

Indications

- Traumatic wounds
- Diabetic foot ulcers
- Pressure ulcers
- Venous ulcers
- Grafts and flaps
- Incisional wounds
- Surgical dehiscence

Contraindications

- Exposed vital organs
- Exposed blood vessels
- Malignant wounds
- Necrotic wounds with eschar
- Untreated osteomyelitis
- Untreated coagulopathy

Evidence shows NPWT promotes faster wound closure and reduces hospital stay compared to Standard Wound Therapy (SWT), leading to quicker recovery and lower healthcare costs.

Precision-Controlled Advanced NPWT by **NoWound**TM



Made in India



CDSO
Licensed



Electrical
Safety



Patent Protected
Design



Cutting Edge Technology

- **Accurate, Regulated Pressure**

Advanced PID Controller: Continuously stabilizes pressure based on wound specifics, eliminating common fluctuations.

Certified Precision: has pressure accuracy certificate from NABL approved lab.

- **Double Tubing**



Maintains optimal exudate removal and pressure, minimizing blockage risks.

- **Real-Time Leakage Monitoring**

Pre-Leakage Indicator: Alerts with green, yellow, and red signals for proactive adjustments, reducing full-scale alarms.

- **Premium dressing kit**

Medical-grade foam & dual-drape dressing kit for optimal coverage, comfort, and accelerated healing.

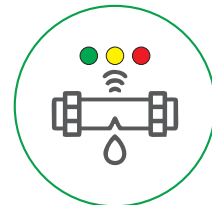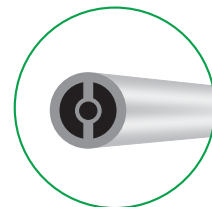

Integrated Care EcoSystem

- **24/7 Helpline**

Expert support available anytime.

- **Automatic Updates**

System-sent therapy progress to clinicians.

- **Professional Assistance**

Daily support from trained staff.

- **Video Guides**

Easy guidance for patients and caregivers.

Designed for Simplicity

- **Portable & Lightweight**

Integrated canister for easy transport.

- **Single-Screen Operation**

Real-time therapy information

- **Car Charger Compatible**

Ideal for mobile patients.

NoWound™ Advanced NPWT

Have a Track Record of Clinical Success

RELIABLE | EFFECTIVE | SEAMLESS | ECONOMICAL

Specifications

200 x 190 x 155 mm size	<1.5 kg weight	12 V DC operating voltage
-40 to -200 mmHg size pressure range	<30 db noise level	8hrs battery backup
Continuous & Intermittent modes		

Accessories



Dressing foam
(Available in three sizes)
S - 7.5x10x3 cm
M - 12.5x18x3 cm
L - 18x25x3 cm



Adhesive
Drape



Scale for
Wound Measurement



Tubing with
Circular Pad



Canister
(600 ml)

24/7 Patient Helpline

+91 9773941323

Read our Clinical Publications



NoWound™
Smart Healing

MEDVITAL VENTURES PVT. LTD.
Second Floor, A-131, Vikas Puri, New Delhi- 110018
GSTIN: 07AARCM8651R1ZB
MD License No: RMD/DCD/24/ZO-5995/1147

Technology Partner:
Inochi Care Private Limited
C-34, Second floor, DDA sheds, Okhla
Industrial Area, Phase 1, Delhi - 110020
License No. MFG/MD/2023/000912

Ph: +91 9773941323 | E: medvitalventures@gmail.com | Web: www.themedvital.com